

## Cartilage confined compression tests

1. Tests include using PBS, longer rest times between tests (sample kept in fridge during rest) and starting preconditioning and stress relaxation 300 microns off of the 10g find contact position to capture data increase from the instant any contact is established.
2. Sample stabilized in PBS for 30 min before thickness measurement and testing.
3. All protocols updated before testing started.
4. Indentor with 3 mm diameter used.
5. Cartilage strip thawed twice.


day 1: oks00TR6-FMC-LAXX-01-01 : thickness: 1.956 mm; Positions: 49.832, 50.0105.

day 2: oks00TR6-FMC-LAXX-01-01 : thickness: 1.928 mm; Positions: 59.0005, 59.2135.


day 3: oks00TR6-FMC-LAXX-01-03 : thickness: 1.88 mm; Positions: 43.7335, 43.9005.

day 4: oks00TR6-FMC-LAXX-01-04 : thickness: 1.914 mm; Positions: 46.7725,46.9855.

day 5: oks00TR6-FMC-LAXX-01-05 : thickness: 1.90 mm; Positions: 55.4788,55.5905.

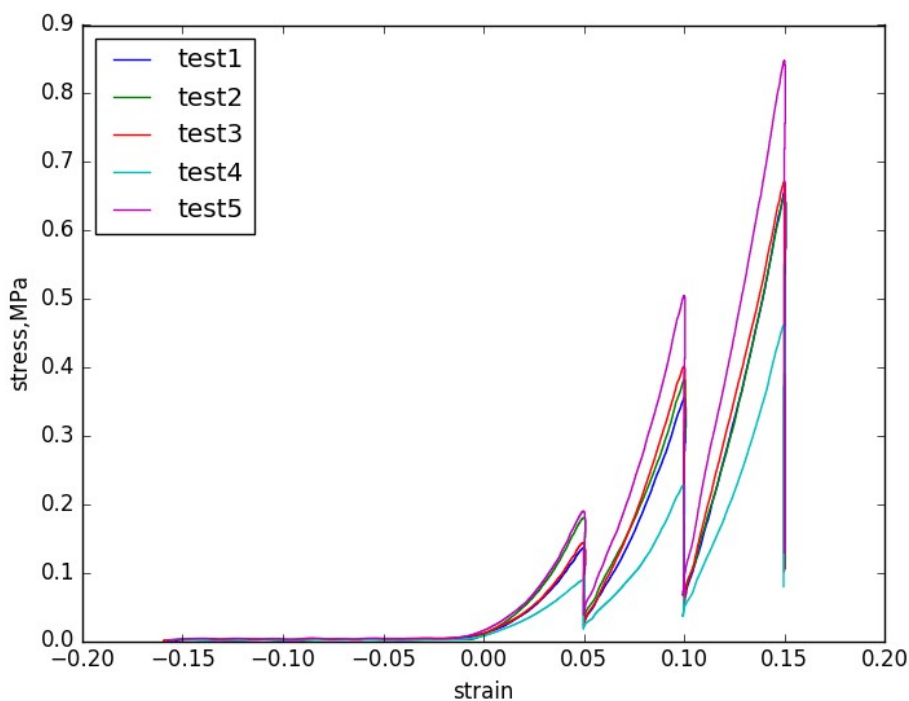

10g force location used as starting position

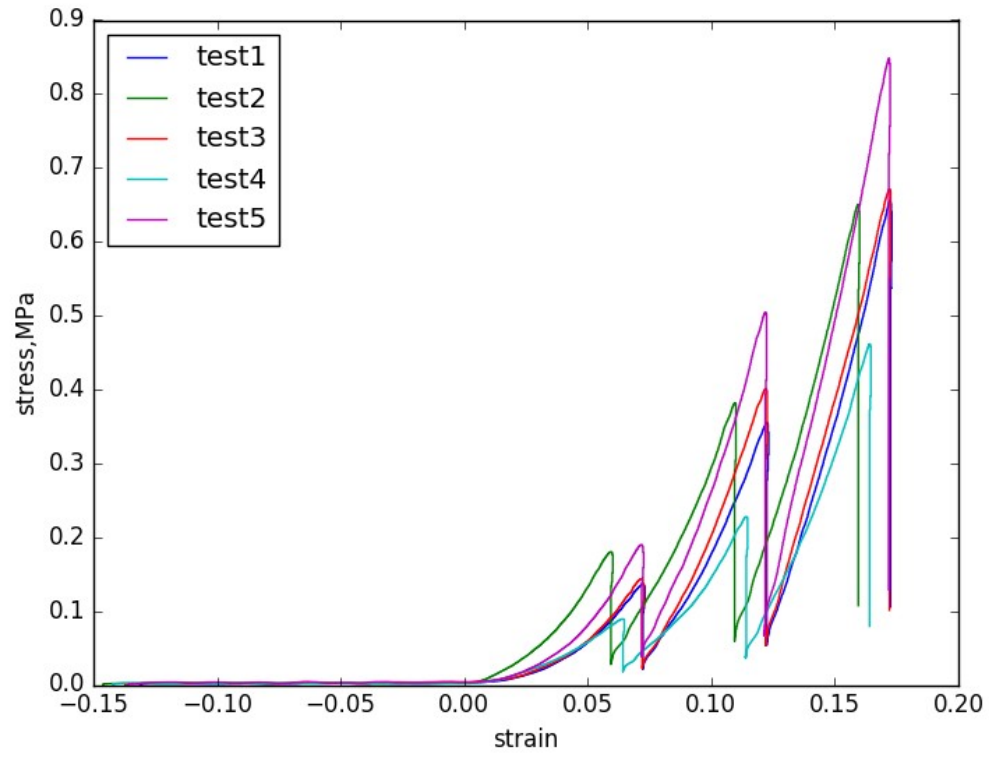

Adjusted to actual zero force-displacement